

INBAUTEK



“Laboratoy table in C made with high resistance compact phenolic resin worktop

Laboratory table in C made with high resistance compact phenolic resin Fundermax Resistance² white/soft grey, 20 mm thick. Ideal for all type of laboratories and indispensable in physical-chemical laboratories. Equipped with metal frame tube cold rolled steel protected with antirust epoxy resins. With a maximum load capacity of 350 kg/linear meter. Four adjustable support points to level the table. Price does not include assembly. Warranty 5 years. Delivery time 4 weeks

reference	dimensions (length x width x height)	pcs/pack	per sales unit
LABT-C01-001	900 x 750 x 900 mm	1	
LABT-C02-001	1200 x 750 x 900 mm	1	
LABT-C03-001	1500 x 750 x 900 mm	1	
LABT-C04-001	1800 x 750 x 900 mm	1	

INBAUTEK



“Laboratory table in A made with high resistance compact phenolic resin worktop

Laboratory table in A made with high resistance compact phenolic resin Fundermax Resistance² white/soft grey, 20 mm thick. Ideal for all type of laboratories and indispensable in physical-chemical laboratories. Equipped with metal frame tube cold rolled steel protected with antirust epoxy resins. With a maximum load capacity of 350 kg/linear meter. Four adjustable support points to level the table. Price does not include assembly. Warranty 5 years. Delivery time 4 weeks

reference	dimensions (length x width x height)	pcs/pack	per sales unit
LABT-A01-001	900 x 750 x 900 mm	1	
LABT-A02-001	1200 x 750 x 900 mm	1	
LABT-A03-001	1500 x 750 x 900 mm	1	
LABT-A04-001	1800 x 750 x 900 mm	1	

INBAUTEK



“Double table for laboratory ovens, incubators and furnaces

Table for furnaces, ovens and incubators at two heights with a separation of 54 cm. Worktop of high resistance compact phenolic resin Fundermax Resistance² white/soft grey, 20 mm thick. Ideal for all type of laboratories. Equipped with metal frame tube cold rolled steel, protected with antirust epoxy resins. Four adjustable support points to level the table. Price does not include assembly. Warranty 5 years. Delivery time: 2 weeks

reference	dimensions (length x width x height)	pcs/pack	per sales unit
LABT-H01-001	900 x 750 x 1200 mm	1	
LABT-H02-001	1200 x 750 x 1200 mm	1	
LABT-H03-001	1500 x 750 x 1200 mm	1	

INBAUTEK



“Auxiliar laboratory cabinet with wheels and one door

Color white. Sides, base and lid made of melamine 19 mm thick. Fronts edged with aluminium 1 mm thick. Right / left door. Hinges with 270° opening. It incorporates a height adjustable shelf, made of melamine 19 mm thick. Four wheels at the base, two with brake. Price does not include assembly. Warranty 5 years. Delivery time: 4 weeks

reference	height of work	dimensions (length x width x height)	pcs/pack	per sales unit
LABC-751-001	750 mm	450 x 520 x 635 mm	1	
LABC-901-001	900 mm	600 x 520 x 805 mm	1	



INBAUTEK

"Auxiliar laboratory cabinet with wheels, one drawer and one door

Color white. Sides, base and lid made of melamine 19 mm thick. Fronts edged with aluminium 1 mm thick. Right / left door. Hinges with 270 ° opening. Drawers with base and bottom made of melamine and guides and sides of sheet steel coated with antioxidant protection. Cushioned closing. Handle made of natural aluminium. It incorporates a height adjustable shelf, made of melamine 19 mm thick. Four wheels at the base, two with brake. Price does not include assembly. Warranty 5 years. Delivery time: 4 weeks

reference	height of work	dimensions (length x width x height)	pcs/pack	per sales unit
LABC-752-001	750 mm	450 x 520 x 635 mm	1	
LABC-902-001	900 mm	600 x 520 x 805 mm	1	



INBAUTEK

"Auxiliar laboratory cabinet with wheels and drawers

Color white. Sides, base and lid made of melamine 19 mm thick. Fronts edged with aluminium 1 mm thick. Drawers with base and bottom made of melamine and guides and sides of sheet steel coated with antioxidant protection. Cushioned closing. Handle made of natural aluminium. It incorporates a height adjustable shelf, made of melamine 19 mm thick. Four wheels at the base, two with brake. Price does not include assembly. Warranty 5 years. Delivery time: 4 weeks

reference	height of work	drawers n°	dimensions (length x width x height)	pcs/pack	per sales unit
LABC-753-001	750 mm	3	450 x 520 x 635 mm	1	
LABC-903-001	900 mm	4	600 x 520 x 805 mm	1	



INBAUTEK

"Auxiliar laboratory cabinet with wheels and two doors

Color white. Sides, base and lid made of melamine 19 mm thick. Fronts edged with aluminium 1 mm thick. Hinges with 270° opening. It incorporates a height adjustable shelf, made of melamine 19 mm thick. Four wheels at the base, two with brake. Price does not include assembly. Warranty 5 years. Delivery time: 4 weeks

reference	height of work	dimensions (length x width x height)	pcs/pack	per sales unit
LABC-904-001	900 mm	900 x 520 x 805 mm	1	




INBAUTEK

"Auxiliar laboratory cabinet with wheels, one drawer and two doors

Color white. Sides, base and lid made of melamine 19 mm thick. Fronts edged with aluminium 1 mm thick. Hinges with 270° opening. Drawers with base and bottom made of melamine and guides and sides of sheet steel coated with antioxidant protection. Cushioned closing. Handle made of natural aluminium. It incorporates a height adjustable shelf, made of melamine 19 mm thick. Four wheels at the base, two with brake. Price does not include assembly. Warranty 5 years. Delivery time: 4 weeks

reference	height of work	dimensions (length x width x height)	pcs/pack	per sales unit
LABC-905-001	900 mm	900 x 520 x 805 mm	1	



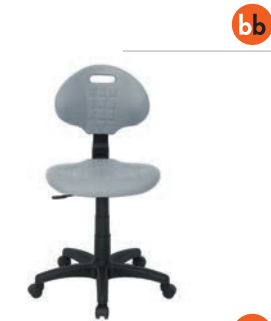


bb

Low laboratory chair with wheels

NEW

Five wheels, rotating, for work in height of 75 cm. Seat/back material: Anti-slip polyurethane, easy to clean. Resistant to abrasion and cleaning products. Adjustment range of seat height: 40-58 cm. Adjustable seat by gas lift. Seat dimensions (width x depth): 46 x 44 cm. Dimensions of the backrest (width x height): 41 x 31 cm. Base: black polyamide with fiberglass, diameter of base: 60 cm. 1 year warranty

reference	description	pcs/pack	per sales unit
LACH-040-001	low laboratory chair with wheels	1	

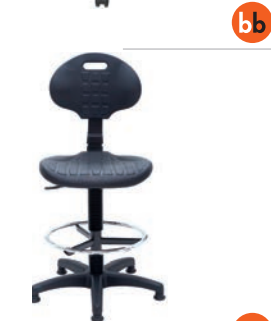


bb

High laboratory chair with footrest

NEW

Five feet, rotating seat, for work in height of 90 cm. Seat/back material: Anti-slip polyurethane, easy to clean. Resistant to abrasion and cleaning products. Adjustment range of seat height: 55-74 cm. Adjustable seat by gas lift. Seat dimensions (width x depth): 46 x 44 cm. Dimensions of the backrest (width x height): 41 x 31 cm. With 45 cm diameter footrest. Base: black polyamide with fiberglass, base diameter: 60 cm. 1 year warranty

reference	description	pcs/pack	per sales unit
LACH-055-001	high laboratory chair with footrest	1	

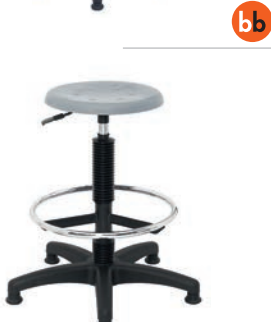

bb

High laboratory stool with footrest

NEW

Five feet, rotating seat. Seat material: Anti-slip polyurethane, easy to clean. Resistant to abrasion and cleaning products. Adjustment range of seat height: 54-85 cm. Adjustable seat by gas lift. Diameter of the seat: 33 cm. With footrest of 45 cm diameter. Base: gray polypropylene, base diameter: 60 cm. 1 year warranty

reference	description	pcs/pack	per sales unit
LASL-054-001	high laboratory stool with footrest	1	

VISIT US AT:

THE WORLD OF SCIENCE AND ANALYTICS @

ARABLAB

DUBAI 2019

26 | 27 | 28

MARS 2019

Paris expo

Porte de Versailles

Hall 4 | FRANCE

MADRID

28-29 March 2019

Hall 8 IFEMA, Stand 146

DUBAI, U.A.E.

12-14 March 2019

PARIS

26-28 March 2019

Hall 4, Stand D44-E45



HEALTH AND SAFETY

Epibox



Visitor safety glass

Made of polycarbonate with anti-scratches treatment and aeration slots in arms. According to EN 166 1.F CE standards

reference	description	pcs/pack	price/pack	average price/pc
SPEC-001-001	visitor safety glass	1		
12 units "Saving package"				
SPEC-001-012	visitor safety glass	12		

Epibox



Safety glass

Made of polycarbonate with anti-scratch treatment. Modern design and ultra light. Model A: According to EN 166 1.F CE standards. Model B: According to EN 166, EN 170 and EN 172 CE standards

reference	description	pcs/pack	price/pack	average price/pc
SPEC-002-001	safety glass (A)	1		
SPEC-UV2-001	safety glass, grey sunscreen with UV filter (B)	1		
12 units "Saving package"				
SPEC-002-012	safety glass (A)	12		

Epibox



Safety goggle

Frame made of comfortable and soft PVC, glass made of polycarbonate. High protection against impacts, particles, splashes, etc. Fixing self-adjustic elastic band. According to EN 166 CE standards

reference	description	pcs/pack	price/pack	average price/pc
SPEC-003-001	safety goggle	1		
10 units "Saving package"				
SPEC-003-010	safety goggle	10		



Safety goggle *Premium Line, LITE series*

Extremely comfortable, resistant frame. Small, ideal for workers who do not need prescription glasses. Forehead section with air vents. Adjustable headband. Optical lens: class 1, antifog treatment and anti-scratch treatment. Medium-energy impact resistance to high-speed particles, protection against large dust particles > 5 microns, protection against liquid droplets. According to EN 166 CE standards

reference	description	pcs/pack	price/pack	
SPEC-004-001	safety goggle, LITE series	1		



Safety reading glasses *Premium Line, WORK&FUN series*

Extremely comfortable, resistant frame. 5 models available. PC lenses with add powers of +1.00 +1.50 +2.00 +2.50 +3.00 diopters and optical lens: class 1. Low-energy impact resistance to high-speed particles. According to EN 166 CE standards. Glasses cord included

reference	diopters	frame color	price/pack	average price/pc
SPEC-D10-001	1	green	1	
SPEC-D15-001	1,5	blue	1	
SPEC-D20-001	2	yellow	1	
SPEC-D25-001	2,5	grey	1	
SPEC-D30-001	3	red	1	



Epibox



Face protector

Face protector made of polycarbonate with head self-adjusting elastic band. Internal foam strip. Very convenient front and side protection (protector covering 180°). High protection against impacts, particles, splashes, etc. According to EN 166 CE standards

reference	description	pcs/pack	price/pack	
FACE-001-001	face protector	1		

Disposable respirator, FFP1 level

Protection against solid particles and low toxicity liquids in low concentration (according to FFP1 level protection defined by EN149-2001 standards); up to 4 x TLV. Spherical shape with adjustable nose clip

reference	type	pcs/pack	price/pack	average price/pc
MASK-001-020	FFP1 spherical shape without valve	20		
MASK-0V1-010	FFP1 spherical shape with valve	10		

Disposable respirator, FFP2 level

Protection against solid particles and low or medium toxicity liquids (according to FFP2 level protection defined by EN149-2001 standards); up to 12 x TLV. Available in spherical shape with valve and folded shape without valve

reference	type	pcs/pack	price/pack	average price/pc
MASK-002-010	FFP2 spherical shape with valve	10		
MASK-P02-020	FFP2 folded without valve	20		

Disposable respirator, with valve, FFP3 level

Protection against solid particles and high toxicity liquids (according to FFP3 level protection defined by EN149-2001 standards); up to 50 x TLV. Folded, with exhalation valve and adjustable nose clip

reference	type	pcs/pack	price/pack	average price/pc
MASK-003-020	FFP3 folded with valve	20		

Disposable cap

Disposable cap made of white polypropylene, circular shape

reference	type	pcs/pack	price/pack	average price/pc
CAPP-001-100	white color	100		

Disposable tie-on mask

Tie-on mask, made of 1 single layer of paper with self-adjusting elastic band

reference	type	pcs/pack	price/pack	average price/pc
MASK-000-100	white color	100		

Disposable visitor coat

Made of polypropylene with pocket on the chest and velcro closure

reference	type	pcs/pack	price/pack	average price/pc
COAP-001-010	one size, white color	10		

Disposable shoe cover

Made of polyethylene

reference	type	pcs/pack	price/pack	average price/pc
SHOP-001-100	green color	100		
SHOP-002-100	blue color, non-slip soles	100		

Epibox



Epibox



Epibox



Epibox



Epibox



Lab-coat for women

With collar flap, closure by 5 buttons and two pockets in the lower side, one on the chest. White color

reference	size	waist girth	chest girth	pcs/pack	price/pack
In 65% polyester/35% cotton					
GOSW-00S-001	S	82 - 88 cm	96 - 102 cm	1	
GOSW-00M-001	M	90 - 96 cm	104 - 110 cm	1	
GOSW-00L-001	L	98 - 104 cm	112 - 118 cm	1	
GOSW-0XL-001	XL	106 - 112 cm	120 - 126 cm	1	
GOSW-XXL-001	XXL	114 - 120 cm	128 - 134 cm	1	
In 65% polyester/35% cotton, "saving package"					
GOSW-00S-010	S	82 - 88 cm	96 - 102 cm	10	
GOSW-00M-010	M	90 - 96 cm	104 - 110 cm	10	
GOSW-00L-010	L	98 - 104 cm	112 - 118 cm	10	
GOSW-0XL-010	XL	106 - 112 cm	120 - 126 cm	10	
GOSW-XXL-010	XXL	114 - 120 cm	128 - 134 cm	10	
In 100% cotton					
GOCW-00S-001	S	82 - 88 cm	96 - 102 cm	1	
GOCW-00M-001	M	90 - 96 cm	104 - 110 cm	1	
GOCW-00L-001	L	98 - 104 cm	112 - 118 cm	1	
GOCW-0XL-001	XL	106 - 112 cm	120 - 126 cm	1	
GOCW-XXL-001	XXL	114 - 120 cm	128 - 134 cm	1	
In 100% cotton, "saving package"					
GOCW-00S-010	S	82 - 88 cm	96 - 102 cm	10	
GOCW-00M-010	M	90 - 96 cm	104 - 110 cm	10	
GOCW-00L-010	L	98 - 104 cm	112 - 118 cm	10	
GOCW-0XL-010	XL	106 - 112 cm	120 - 126 cm	10	
GOCW-XXL-010	XXL	114 - 120 cm	128 - 134 cm	10	

Epibox



Lab-coat for men

With collar flap, closure by 5 buttons and two pockets in the lower side and one on the chest. White color

reference	size	chest girth	pcs/pack	price/pack
In 65% polyester/35% cotton				
GOSM-00S-001	S	98 - 104 cm	1	
GOSM-00M-001	M	106 - 112 cm	1	
GOSM-00L-001	L	114 - 120 cm	1	
GOSM-0XL-001	XL	122 - 128 cm	1	
GOSM-XXL-001	XXL	130 - 136 cm	1	
In 65% polyester/35% cotton, "saving package"				
GOSM-00S-010	S	98 - 104 cm	10	
GOSM-00M-010	M	106 - 112 cm	10	
GOSM-00L-010	L	114 - 120 cm	10	
GOSM-0XL-010	XL	122 - 128 cm	10	
GOSM-XXL-010	XXL	130 - 136 cm	10	
In 100% cotton				
GOCM-00S-001	S	98 - 104 cm	1	
GOCM-00M-001	M	106 - 112 cm	1	
GOCM-00L-001	L	114 - 120 cm	1	
GOCM-0XL-001	XL	122 - 128 cm	1	
GOCM-XXL-001	XXL	130 - 136 cm	1	
In 100% cotton, "saving package"				
GOCM-00S-010	S	98 - 104 cm	10	
GOCM-00M-010	M	106 - 112 cm	10	
GOCM-00L-010	L	114 - 120 cm	10	
GOCM-0XL-010	XL	122 - 128 cm	10	
GOCM-XXL-010	XXL	130 - 136 cm	10	



Laboratory scrub set

Set of unisex laboratory scrub with long sleeve top. Made of twill (65% polyester/35% cotton), white. Includes a V-neck top with a chest pocket and two waist patch pockets, as well as a pair of elastic waistband pants with 2 italian pockets. Grammage: 195 g

reference	size	waist girth	pcs/pack	price/pack
LSSS-OXS-001	XS	72 cm	1	
LSSS-OOS-001	S	76 - 80 cm	1	
LSSS-OOM-001	M	84 - 88 cm	1	
LSSS-OOL-001	L	92 - 96 cm	1	
LSSS-OXL-001	XL	100 - 104 cm	1	
LSSS-XXL-001	XXL	108 - 112 cm	1	
LSSS-3XL-001	XXXL	116 - 120 cm	1	
10 units "Saving package"				
LSSS-OXS-010	XS	72 cm	10	
LSSS-OOS-010	S	76 - 80 cm	10	
LSSS-OOM-010	M	84 - 88 cm	10	
LSSS-OOL-010	L	92 - 96 cm	10	
LSSS-OXL-010	XL	100 - 104 cm	10	
LSSS-XXL-010	XXL	108 - 112 cm	10	
LSSS-3XL-010	XXXL	116 - 120 cm	10	

Nitrile disposable gloves

Powder-free, non sterile. Nitrile is much more chemical resistant than other kind of disposable gloves commonly used in laboratory, such as latex or vinyl. 100% latex protein free: An alternative to common latex gloves to prevent allergies. Ambidextrous. AQL 1,5. According to 89/686/CEE, EN 420, EN 374 1/2/3 and EN 388 standards

reference	type	pcs/pack	price/pack
Powder-free			
GLON-OOS-100	powder free, S (6-7) size	100	
GLON-OOM-100	powder free, M (7-8) size	100	
GLON-OOL-100	powder free, L (8-9) size	100	
GLON-OXL-100	powder free, XL (9-10) size	100	
10 x 100 Powder-free gloves "Saving package"			
GLON-OOS-1K0	powder free, S (6-7) size	10 x 100	
GLON-OOM-1K0	powder free, M (7-8) size	10 x 100	
GLON-OOL-1K0	powder free, L (8-9) size	10 x 100	
GLON-OXL-1K0	powder free, XL (9-10) size	10 x 100	

Nitrile disposable gloves "Sensitive"

Powder-free, non sterile. More sensitive to touch than usual nitrile gloves. With the advantages of nitrile gloves with a more competitive price. 100% latex protein free. Ambidextrous. AQL 1,5. According to 89/686/CEE, EN 420, EN 374 1/2/3 and EN 388 standards

reference	type	pcs/pack	price/pack
Powder-free			
GLOS-OOS-100	powder free, S (6-7) size	100	
GLOS-OOM-100	powder free, M (7-8) size	100	
GLOS-OOL-100	powder free, L (8-9) size	100	
GLOS-OXL-100	powder free, XL (9-10) size	100	
10 x 100 Powder-free gloves "Saving package"			
GLOS-OOS-1K0	powder free, S (6-7) size	10 x 100	
GLOS-OOM-1K0	powder free, M (7-8) size	10 x 100	
GLOS-OOL-1K0	powder free, L (8-9) size	10 x 100	
GLOS-OXL-1K0	powder free, XL (9-10) size	10 x 100	

Epibox

Epibox

Epibox


Epibox



Latex disposable exam gloves

Available in powdered and powder-free, non sterile. 100% Pure natural latex, natural colour. Ambidextrous. AQL 1.5. According to 89/686/CEE, EN 420 and EN 374 1/2 standards

reference	type	pcs/pack	price/pack
Powdered			
GLOL-PXS-100	powdered, XS (5-6) size	100	
GLOL-P0S-100	powdered, S (6-7) size	100	
GLOL-P0M-100	powdered, M (7-8) size	100	
GLOL-P0L-100	powdered, L (8-9) size	100	
GLOL-PXL-100	powdered, XL (9-10) size	100	
Powder-free			
GLOL-0XS-100	powder free, XS (5-6) size	100	
GLOL-00S-100	powder free, S (6-7) size	100	
GLOL-00M-100	powder free, M (7-8) size	100	
GLOL-00L-100	powder free, L (8-9) size	100	
GLOL-0XL-100	powder free, XL (9-10) size	100	
10 x 100 Powdered gloves "Saving package"			
GLOL-PXS-1K0	powdered, XS (5-6) size	10 x 100	
GLOL-P0S-1K0	powdered, S (6-7) size	10 x 100	
GLOL-P0M-1K0	powdered, M (7-8) size	10 x 100	
GLOL-P0L-1K0	powdered, L (8-9) size	10 x 100	
GLOL-PXL-1K0	powdered, XL (9-10) size	10 x 100	
10 x 100 Powder-free gloves "Saving package"			
GLOL-0XS-1K0	powder free, XS (5-6) size	10 x 100	
GLOL-00S-1K0	powder free, S (6-7) size	10 x 100	
GLOL-00M-1K0	powder free, M (7-8) size	10 x 100	
GLOL-00L-1K0	powder free, L (8-9) size	10 x 100	
GLOL-0XL-1K0	powder free, XL (9-10) size	10 x 100	

Epibox



Autoclave gloves

Made of Nomex III fibers (mixture of thermo-resistant and infusible aromatic polyamides) and with no asbestos fibers. With normal service conditions, physical properties stay unchanged under several hours at 260 °C, withstanding as well 370 °C when its degradation starts. When approached to flame, they do not melt, but instead they get scorched making a crust allowing protection. Can be used as well to handle sporadically items stored under nitrogen vapour phase (-180 °C). Resistant to bases, acids and organic products but not to large exposures to some concentrated inorganics acids.

Available in 30 cm long or 52 cm long. According to EN 407 y EN 388 standards

reference	type	pair of gloves	price/pack
GLOT-S30-001	glove shape, small size, 30 cm long	1	
GLOT-M30-001	glove shape, medium size, 30 cm long	1	
GLOT-L30-001	glove shape, large size, 30 cm long	1	
GLOT-S52-001	glove shape, small size, 52 cm long	1	
GLOT-M52-001	glove shape, medium size, 52 cm long	1	
GLOT-L52-001	glove shape, large size, 52 cm long	1	

Epibox



Anti-acid PVC gloves

Gloves made of red PVC with smooth finishing. Very good resistance to abrasion, concentrated acids, hydrocarbons, oils and grease. According to EN 388 standards

reference	type	pair of gloves	price/pack
GLOP-035-001	one size, length 40 cm	1	

Cryogenic gloves

Used to handle items stored under liquid nitrogen or other very cold temperatures (from -120 °C up to -190 °C). Designed in order to leave plenty of space to fingers. Water-resistant. Up to mid-arm, 37 cm long. According to EN 511 and EN 388 standards

reference	type	pair of gloves	price/pack
GLOF-S37-001	up to mid-arm, small size	1	
GLOF-M37-001	up to mid-arm, medium size	1	
GLOF-L37-001	up to mid-arm, large size	1	

Earmuff

Very light earmuff. Light noise reduction (SNR = 27 db). Adjustable height pads and padded harness. According to EN 352-1 CE standards

reference	type	pcs/pack	price/pack
EARP-001-001	SNR = 27 dB	1	

Safe-Handle bottle protector

Prevents breakage of 2.5 L glass bottles commonly used for most laboratory reagents and solvents manufacturers. If a bottle does break, it helps contain glass shards and liquid spills, protecting lab personnel. Integrated handle helps in handling and pouring heavy bottles. Easy-to-use mechanism and fully reusable. Important: suitable only for bottle models with smooth external surface (not suitable for bottles with integrated handle)

reference	description	pcs/pack	price/pack
SAFE-001-001	Safe-Handle bottle protector	1	

Epibox



Epibox



Epibox



plum



Eyewash stationary stand

Wall stand with two eyewash bottles, explanatory pictograms for instructions and a small mirror. Specially designed for work areas with eye damage risks by chemical products or foreign bodies. Includes two bottles available separately for sale: one including 500 ml of eyewash liquid (sterile NaCl 0.9%) and one including 200 ml of neutral solution (sterile phosphate buffer 4.9%). Neutral solution is to be used to neutralize quickly acids and bases and must be rinsed with the eyewash solution afterwards. Both bottles have an ergonomic mouth perfectly adapted to the eye, according to CE standards and with expiry date from 2 to 3 years

reference	description	pcs/pack	price/pack
EYES-S01-001	stationary stand eyewash	1	
Spare parts			
EYES-SAL-500	saline solution, 500 ml	1	
EYES-NEU-200	200 ml neutral pH solution	1	

plum



DUO Eyewash stationary stand

Wall stand with two DUO eyewash bottles, explanatory pictograms for instructions and a small mirror. Specially designed for work areas with eye damage risks by chemical products or foreign bodies. Includes two bottles available separately for sale: one including 500 ml of eyewash liquid (sterile NaCl 0.9%) and one including 200 ml of neutral solution (sterile phosphate buffer 4.9%). Neutral solution is to be used to neutralize quickly acids and bases and must be rinsed with the eyewash solution afterwards. The DUO bottles are equipped with a special eye device that allows rinsing of both eyes at the same time, according to CE standards and with expiry date from 2 to 3 years

reference	description	pcs/pack	price/pack
EYED-S01-001	DUO eyewash stationary stand	1	
Spare parts			
EYED-SAL-1K0	saline solution, 1000 ml, DUO bottle	1	
EYED-NEU-500	500 ml neutral pH solution, DUO bottle	1	



NEW

labkem



Absorbent powder for spilled liquids AUX - Auxiliary grade

Naturally occurring absorbent particles specially designed to absorb all kinds of spills or leaks that occur inside or outside in protected areas of the wind. Mineral sepiolite 30/60 assay > 65%, max. bulk density 720 g/l, min. water absorption 110% and min. oil absorption 65%.
INSTRUCTIONS: Add the granules on the spilled liquid, wait a few minutes until the liquid is completely absorbed and collect it with a shovel. For further information, see Specification Sheet at www.labbox.com




reference	package	price/pack
ABSU-00P-1K0	 1 kg	
ABSU-00P-5K0	 5 kg	

labkem



Detergents for laboratory use, for manual cleaning

Recommended for manual use, soaking and ultrasonic cleaning of glassware and laboratory equipment made of stainless steel, porcelain, ceramic and plastic. LABKEM CLEANER is biodegradable and ensures complete removal: grease, oils, blood, cell tissue, chemicals and biological residues, etc. Complete rinse without residues. See Specification Sheet to consult the composition and other information (to be downloaded at www.labbox.com)

reference	description	package	price/pack
SOAP-066-5K0	Labkem cleaner M66 Alkaline liquid detergent	 5 L	
SOAP-067-5K0	Labkem cleaner M67 Alkaline liquid detergent, Phosphate free	 5 L	
SOAP-068-5K0	Labkem cleaner M68 Neutral liquid detergent	 5 L	

Detergents for laboratory use, for automatic washing machine

labkem



reference	description	package	price/pack
SOAP-102-5K0	Labkem cleaner A102 Alkaline liq. detergent, Phosphate free	5 L	
SOAP-106-5K0	Labkem cleaner A106 Neutral liquid detergent, Phosphate free	5 L	
SOAP-107-5K0	Labkem cleaner A107 Alkaline powder detergent	4,5 kg	

Rinsing agent, for automatic washing machine

labkem

Acidic neutralizing agent - surfactant free for use in laboratory glassware washers. LABKEM RINSE AGENT A103 is biodegradable and is used in the neutralizing cycle after the main alkaline washing cycle with LABKEM CLEANER A102. It neutralizes and removes residual alkaline and mineral traces which may cause white stains on glassware. LABKEM RINSE AGENT A103 can also be used in the prewash cycle. Complete rinse without residues. See Specification Sheet to consult the composition and other relevant information (to be downloaded at www.labbox.com)



reference	description	package	price/pack
SOAP-103-5K0	labkem rinse agent A103	5 L	

Disinfectants for laboratories

labkem

Disinfectants for surfaces, material and laboratory equipment. 3 varieties available for different types of laboratories: General, Medical and Molecular Biology. **DSNF-GFS-750** - GERMfree solution for disinfection of instruments and surfaces in general. **DSNF-MYS-250** - MYCOPLASMAfree solution for particularly active disinfection of mycoplasma and a wide variety of microorganisms. **DSNF-ZDS-250** - ZERO DNA solution for cleaning contaminating DNA and RNA particles of examinations, samples and tests in PCR assays. Complete rinse without residues. See Specification Sheet to consult the composition and other relevant information (to be downloaded at www.labbox.com)



reference	description	package	price/pack
DSNF-GFS-750	GERMfree solution	750 ml	
DSNF-GFR-750	GERMfree refill	750 ml	
DSNF-MYS-250	MYCOPLASMAfree solution spray	250 ml	
DSNF-MYR-500	MYCOPLASMAfree refill	500 ml	
DSNF-ZDS-250	ZERO DNA solution spray	250 ml	
DSNF-ZDR-500	ZERO DNA refill	500 ml	

Draining rack

bb

PVC coated wire. Can be hanged or placed upon a table. Includes a water collection tray




reference	type	h (mm)	w (mm)	base (mm)	pcs/pack	price/pack
DRAI-055-001	55 pos.	650	360	360 x 130	1	

NEW

Wall draining rack for labware

Made of ABS and polycarbonate. To be used hanging on the wall. Bottom of the rack is provided with a tray which collects all the waste water with drain outlet to facilitate emptying

bb

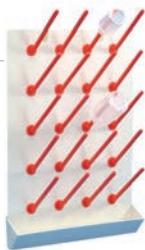

reference	type	h (mm)	w (mm)	pcs/pack	price/pack
DRAI-P20-001	20 pos	460	350	1	

bb



astik's

Test tube basket

Made of polypropylene. Autoclavable at 121°C

reference	h (mm)	l (mm)	w (mm)	pcs/pack	price/pack
TTBP-001-001	140	120	110	1	
TTBP-002-001	160	160	160	1	

Waste disposal bag resistant to autoclave, HDPE

Made of polyethylene (HDPE). Transparent. Perfect for the sterilization and disposal of hazardous and contaminated labware. Thickness: 0,03 mm

reference	capacity	dimensions (mm)	pcs/pack	price/pack	average price/pc
WABS-300-100	8 L	300 x 500	100		
WABS-400-100	38 L	400 x 780	100		
WABS-600-100	60 L	600 x 800	100		
"Saving package"					
WABS-300-2K0	8 L	300 x 500	20 x 100		
WABS-400-1K0	38 L	400 x 780	10 x 100		
WABS-600-500	60 L	600 x 800	5 x 100		



astik's

Waste disposal bag resistant to autoclave, PP

Made of polypropylene (PP), autoclavable at 134°C. Perfect for the sterilization and disposal of hazardous and contaminated labware. For autoclaving, the bag must not be completely closed. Thickness: 0,05 mm

reference	capacity	dimensions (mm)	pcs/pack	price/pack	average price/pc
WABP-300-100	8 L	300 x 500	100		
WABP-400-100	38 L	400 x 780	100		
WABP-600-100	60 L	600 x 800	100		
"Saving package"					
WABP-300-1K5	8 L	300 x 500	15 x 100		
WABP-400-600	38 L	400 x 780	6 x 100		
WABP-600-300	60 L	600 x 800	3 x 100		



bb

Stainless steel wire basket

Designed for transporting, storing, cleaning and sterilising laboratory utensils. Excellent chemical resistance. Both models round and square comes with folding handle, and the square model comes also with lid

reference	external dimensions (length x width x height)	type	pcs/pack	price/pack
BASS-R15-001	Ø150 x h120 mm	round	1	
BASS-R21-001	Ø210 x h180 mm	round	1	
BASS-R24-001	Ø240 x h180 mm	round	1	
BASS-R30-001	Ø300 x h200 mm	round	1	
BASS-S16-001	160 x 160 x 160 mm	square	1	
BASS-S20-001	200 x 200 x 200 mm	square	1	
BASS-S30-001	300 x 200 x 150 mm	square	1	



bb

Burette brushes

Natural bristle hand on twisted galvanised iron wire stem

reference	l _{total} (mm)	Ø _{brush} (mm)	l _{brush} (mm)	type	pcs/pack	price/pack
BRU3-005-001	800	25	150	natural bristle	1	



NEW

Test tube brushes

White nylon or natural bristle hand on twisted galvanised iron wire stem

reference	l _{total} (mm)	Ø _{brush} (mm)	l _{brush} (mm)	type	pcs/pack	price/pack	average price/pc
BRU1-002-002	120	6	40	nylon (A)	2		
BRU3-001-001	170	17	55	natural bristle (B)	1		
BRU3-002-001	200	25	75	natural bristle (C)	1		
BRU3-003-001	240	35	70	natural bristle (D)	1		

Brushes for general purposes

Natural bristle hand on twisted galvanised iron wire stem

reference	l _{total} (mm)	Ø _{brush} (mm)	l _{brush} (mm)	type	pcs/pack	price/pack
BRU4-001-001	280	40	80	natural bristle (A)	1	
BRU4-002-001	310	50	100	natural bristle (B)	1	
BRU4-003-001	355	50	110	natural bristle (C)	1	
BRU4-004-001	600	70	150	natural bristle (D)	1	

Brush for special purposes

Natural bristle

reference	l _{total} (mm)	Ø _{brush} (mm)	l _{brush} (mm)	type	pcs/pack	price/pack
BRU2-001-001	380	16	120	natural bristle, for big beakers (A)	1	
BRU2-003-001	580	40	220	natural bristle, for round bottom flasks (B)	1	

Set of 9 assorted brushes

Includes brushes for general use and test tubes of different sizes in nylon and natural bristle

reference	type	pcs/pack	price/pack
BRU1-KIT-001	9 assorted brushes	1	

Industrial work cellulose paper

Double layer. Embossing and precut. Ideal for basic drying tasks. 475 m (1610 sheets)

reference	Ø _{hole}	wide	weight	pcs/pack	price/pack	average price/pc
IWCP-001-002	7 cm	26 cm	4,5 Kg	2		

Dry hands cellulose paper, multipurpose

Double laminated layer. Embossing and precut. Ideal for basic drying tasks. 130 m (440 sheets)

reference	Ø _{hole}	wide	long	weight	pcs/pack	price/pack	average price/pc
DHCP-001-006	7 cm	20 cm	130 m	800 g	6		


NEW

Cellulose Tissues


Double layer. 100% cellulose. 100 service dispenser box. Suitable for all types of environments

reference	l (mm)	w (mm)	surface weight	pcs/pack	price/pack
CELT-001-001	210	200	32 g	1	

bb

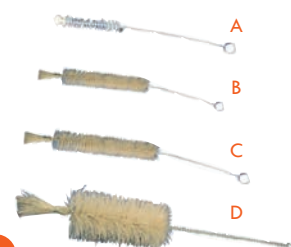

bb

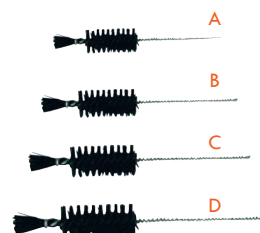

bb



bb



bb



bb



bb



Range of products

BASICHEM Chemicals in general use in laboratories. It includes the following grades:

- AGR** – **Analytical Grade**: for analytical applications with low impurity level (on the order of ppm), according ACS or ISO standards
- EPR** – **Extra Pure**: for ORGANIC SYNTHESIS, CHEMICAL PREPARATIVE or ROUTINE ANALYSIS applications.
- GLR** – **General Laboratory Reagent**: for general use.
- VSA** – **Volumetric Solution for Analysis**: used in volumetric analysis, mostly ACID-BASE.
- IND** – **Indicator**: INDICATORS for volumetric titrations.
- AUX** – **Auxiliary Grade**: auxiliary products for laboratory uses.

LIFEKEM Products for biological and microbiological applications. It includes the following grades:

- MIC** – **Microscopy Grade**: DYES and STAINS in powder or ready to use in microscopy and histology applications.
- HIS** – **Histology Grade**: to use in the field of pathology.
- BAC** – **Bacteriology Grade**: culture media, bases, broths and ingredients for bacteriological applications.
- GEN** – **Molecular Genetic Grade**: culture media, bases and agaroses for MOLECULAR BIOLOGY applications.

CERTIKEM Standard solutions for instrument calibration. It includes a single grade:

- STD** – **Standard grade**: products for calibration of instruments. It includes pH and CONDUCTIVITY

INSTRUKEM High performance products for GC, HPLC, LC-MS and ICP. It includes the following grades:

- GGR** – **Gradient Grade**: high performance solvents for HPLC (High Pressure Liquid Chromatography).
- PGR** – **Pesticide Grade**: solvents for PESTICIDE ANALYSIS.
- LGR** – **LC/MS Grade**: ultra-high performance solvents for LC-MS (Liquid Chromatography – Mass Spectrometry).
- TGR** – **Trace Grade**: very high purity ACIDS, with an impurity level on the order of 0,1-1 ppb, for sample preparation for TRACE ANALYSIS by ICP-MS/OES (Inductively Coupled Plasma – Mass Spectrometry / Optical Emission Spectrometry).

Packing



plastic bottle GL45

plastic bottle

glass bottle GL45

glass bottle

plastic dropper bottle

plastic ampoule

aluminium container



plastic twin neck

Bag in Box PET/Al

glass bottle PP28

metal can

glass dropper bottle

specimen Container

plastic jerrycan



COEX plastic bottle

plastic Bag

COEX plastic jerrycan

plastic pail

vial

bag

zipper bag